



## F. Load Samples

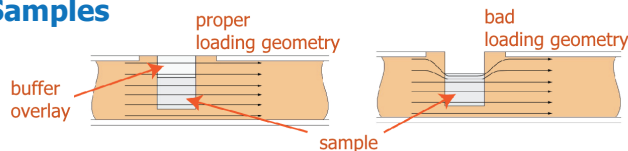

1. **Check the buffer level in the sample wells.** Make sure that the sample wells are completely full to the top with electrophoresis buffer. Fill with additional buffer if necessary.
2. **Remove 40µl of buffer from the first sample well and load 40µl of the sample (or external marker, if assigned to that lane) into that well.** Take care not to pierce the agarose with the pipette tip. There is gel on all sides and bottom of the sample well. In addition, there is an agarose "chimney" surrounding the top of the sample well that protrudes up through the cassette cover. When removing buffer, some users find it useful to immerse the pipette tip just below the surface of the buffer and follow the liquid level down with the tip as the buffer is removed. When adding sample, place the tip of the pipette just below the surface of the buffer, and follow the liquid level up with the tip as the well fills. Don't be concerned if the sample well slightly overfills. The density of the sample will allow it to sink before it can flow out of the well.
3. **Repeat step 2 for the remaining four wells. Make sure the external marker has been loaded into the assigned reference lane.**

## G. Run

1. Close the lid, go to the Main Tab, and make sure the proper protocol is loaded in the "Protocol Name" field.
2. Press "START". The run will automatically stop when every collection is complete.

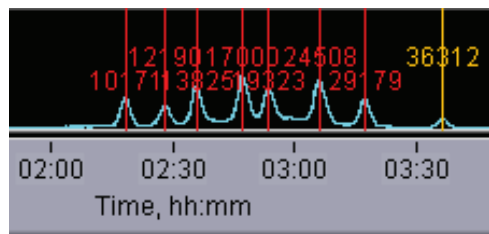

External marker T1 should appear as shown

Approximate Run Times for  
"Tight" Range Targets (hr:min)

kb	hh:mm
18	03:45
23	04:00
27	04:15

## H. Collect Fractions

1. **Remove samples using a standard 100-200µl pipette.** Samples will be in a Tris-TAPS buffer at a volume of 40µl.
2. **Remove the cassette and dispose of properly.** Do not keep used cassettes in the BluePippin with the cover closed. Humidity from the cassette may cause damage to the electrodes.

Help: [support@sagescience.com](mailto:support@sagescience.com) or call 978.922.1832

# BluePippin™ Quick Guide

## 0.75% Agarose Gel Cassette

for targets between 18-27kb

External Marker Only

Product Number: BMF7510

Cassette Definitions: 0.75% DF 18-27kb Marker T1

Cassette and Marker: 0.75% Agarose / Marker T1

Research Use Only

### Recommended Sample Load Guidelines

**Ionic strength:** The ionic strength of the sample should be lower than the ionic strength of the buffer (80mM monovalent ions).

**Protein in the sample:** For best results, samples should be de-proteinized prior to loading if possible.

Maximum Load: 5µg sheared genomic DNA


Minimum Load: low single nanograms

## A. Prepare DNA Samples and External Marker

1. Bring DNA samples up to 30µl with TE.
2. Bring the external marker (**T1**) and the loading solution to room temperature.
3. Combine DNA samples with 10µl of the loading solution.
4. Mix DNA samples thoroughly (vortex mixer). Briefly centrifuge to collect.

## B. Program a Protocol

1. In the BluePippin Software, go to the Protocol Editor Tab.
2. Click "Cassette" folder, and select the appropriate cassette definition.
3. Select the collection mode for each lane (usually "Tight" or "Range"), and enter the size selection parameters.
4. Determine the lane into which the DNA marker (**T1**) will be added and enter it into the "Reference Lane" field. In the example below, **T1** is loaded in lane 1.

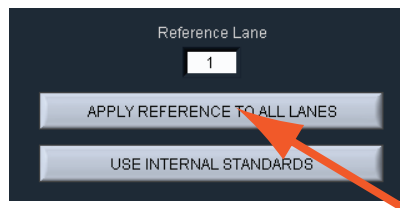

5. Make sure each "Ref Lane" field contains the marker lane designation:

	Tight	Range	Time	Peak	Ref Lane
5					1
4					1
3					1
2					1
1					1

external marker in lane 1

6. Press "Save As" and name and save the protocol.

## C. Calibrate the Optics with the Calibration Fixture

1. Place the calibration fixture onto the optical nest. Close the Lid.
2. Press "CALIBRATE" to launch the calibration window.
3. Press "CALIBRATE" button in the window, and when complete, press "EXIT".




Optical Calibration Fixture

## D. Inspect the Gel Cassette

1. Remove the cassette from the foil packaging.
2. **Inspect the levels of buffer in all buffer reservoirs.** Reservoirs should be nearly full. If the buffer level in any reservoir appears less than 50% full, fill with spare buffer.

anode buffer chambers, side view

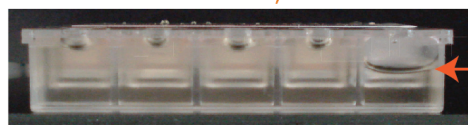

buffer reservoir on right is low fill with spare buffer


3. **Inspect the gel columns.** Look for obvious breakage of the agarose column in each channel. If there is obvious breakage, do not use the lane. Remaining lanes can be used.

bottom of cassette

4. **Inspect for bubbles** due to separation of agarose from bottom of cassette in the region used for optical detection of DNA. If a bubble is observed, do not use the lane. A bubble on the top of the gel column will not affect detection.



optical region

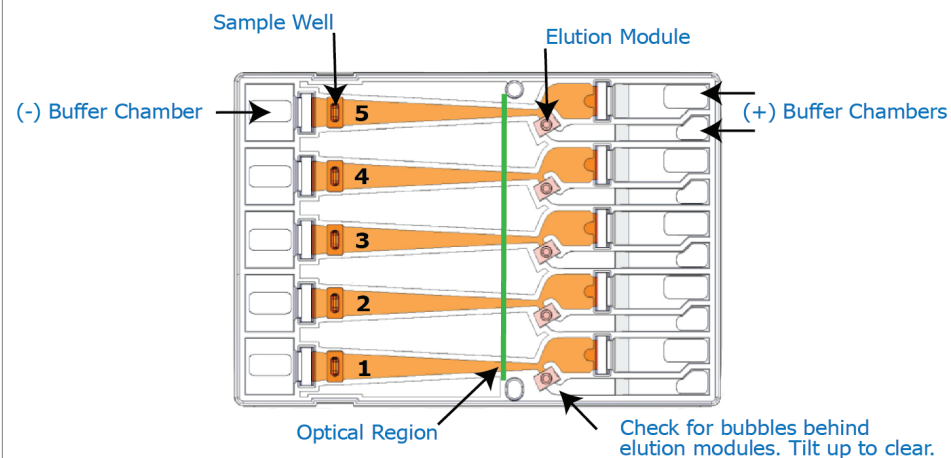

## E. Prepare the Cassette for Loading

1. **Dislodge bubbles from behind the elution modules.** Tilt the cassette, sample well side down, to release any trapped bubbles behind the elution modules.
2. **Place cassette into the optical nest.** Keep the cassette slightly tilted down so that the bubbles in the elution reservoirs don't return to the area behind the elution modules. Be sure the cassette is fully seated into the bottom of the nest to ensure proper optical alignment.
3. **Remove the white tabbed adhesive strips from the cassette.** Place one hand on the cassette and hold it firmly in the nest. Grab the white tabs of the tape and pull the strips firmly and slowly toward the front of the BluePippin until they are removed.
4. **Remove buffer from all elution modules and replace with 40µl of fresh electrophoresis buffer.** Make sure that the pipette tips extend all the way to the bottom of the elution modules without sealing the elution port opening. Test the tip fit using the empty rinse cassette supplied with the instrument.
5. **Seal the elution modules with the adhesive tape strips.** Tape for sealing the elution modules are supplied with cassette packaging. Place tape over the elution modules and rub firmly to fix the tape in position.
6. **Check the buffer level in the sample wells.** Immediately prior to loading, sample wells should be completely filled to the top with buffer. If any wells are underfilled, top them up with additional buffer.
7. **Perform the continuity test.** Close the lid and press the "Test" button located in the lower right area of the Main screen. The test routine runs automatically and measures the current in each separation and elution channel and should return a "PASS" for each separation and elution channel. The cassette temperature must be above 17°C (62°F).